i Neck

3 Pro

breo

breo

CONTENTS

Precautions	03
Product configuration	08
Product features	09
Product functions	10
Battery removal steps	23
Troubleshooting	24
Protection	25
Product information	30

EN

Precautions

Please read this instruction manual carefully before using. The company assumes no responsibility for any consequences caused by misuse or failure of following these instructions.

	Do not use this product while driving.
	Please remove any accessories you might be wearing or other objects that might affect the massaging process before use.
	Do not use the product in the bathroom or other humid conditions.
	Keep the product away from water, fire, and corrosive environment.
	Please comply with relevant safety regulations while using the product under specific circumstances.
	Do not use the product on surgical areas.
	Do not use the product on areas of bruises, damaged skin, swelling, or inflammation.

	Do not use the product if you have cardiovascular (or once had), acute, or dermatologic diseases, malignant tumors, or other abnormal symptoms. Do not use the product in case of hematologic disorders with hemorrhagic symptoms, osteoporosis, soft tissue diseases, or strokes. In case of restricted physical capabilities, physical or mental problems, or mental disabilities, do not use the device without medical consultation first. Stop the massage if you find it unpleasant.
	Do not use the device during menstruation or pregnancy.
	This device is not intended for use by children under the age of 8. Children aged between 8 and 12 must be supervised while using the device.
	Clean with alcohol wipes. Avoid contact with water.

	<p>Avoid contact with sharp objects.</p>
	<p>The device is for private use only and cannot be used as medical equipment. It is not intended for medical diagnosis, treatment, or commercial purposes.</p> <p>Please use only the supplied charging device to charge this product and ensure electrical safety.</p> <p>If any fault or malfunction should occur during use, cut off the power immediately. Never attempt to disassemble or repair the device yourself.</p> <p>Some specific models come with pyretic moxibustion function. The surface of the device heats up during use. Persons insensitive to heat should be especially cautious.</p>

The instructions for class III appliances shall state that it must only be supplied at safety extra low voltage corresponding to the marking on the appliance. This instruction is not necessary for battery-operated appliances if the battery is a primary battery or secondary battery charged outside of the appliance.

EN

NOTES

Necessary knowledge

※ Frequency

1 time per day, no more than 10 minutes.

Overuse might cause discomfort or neck and skin injuries.

When you use the massager for the first time, be mindful of how you feel and observe your shoulders, neck, and skin; do not massage for too long.

※ Size

This product fits persons with neck circumference between 30 and 40 cm.

※ How to use

During use, please close your eyes, lean back, and relax.

Please take off your necklaces, neckhung headphones, ID badges, etc.


Stop the massage and seek medical help immediately if you find it is painful or unpleasant.

※ **After Use**


The skin around your neck may become a little red after massaging. This phenomenon is normal and will quickly disappear.


※ **Others**

The massage head might buzz when it rotates and changes position. This is not a sign of defect or malfunction.

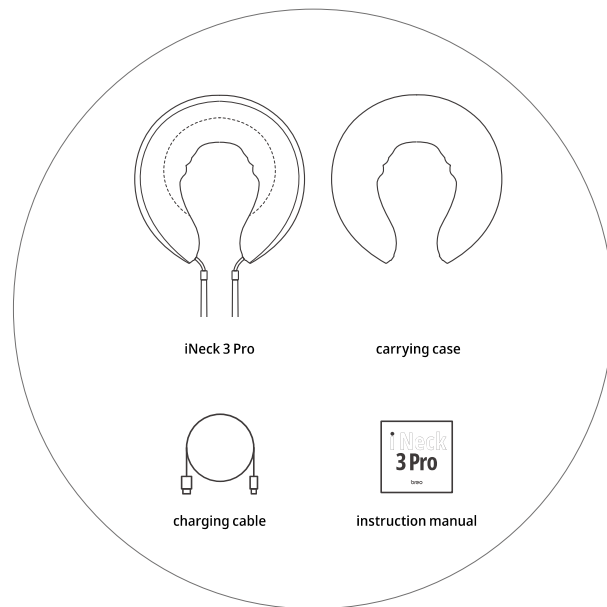PRODUCT CONFIGURATION

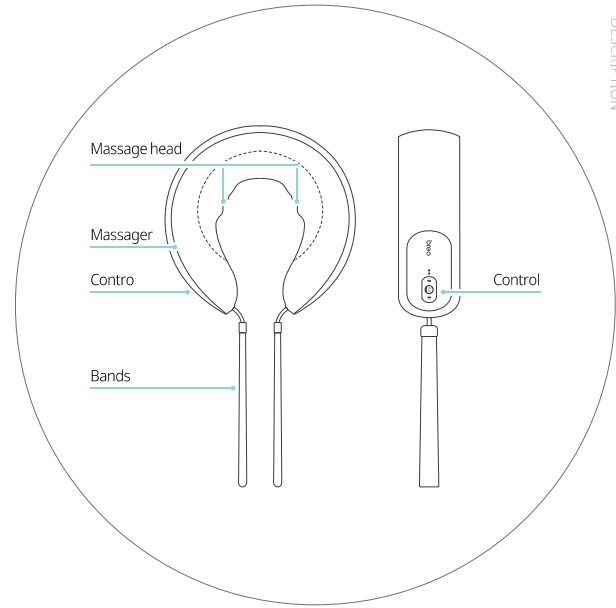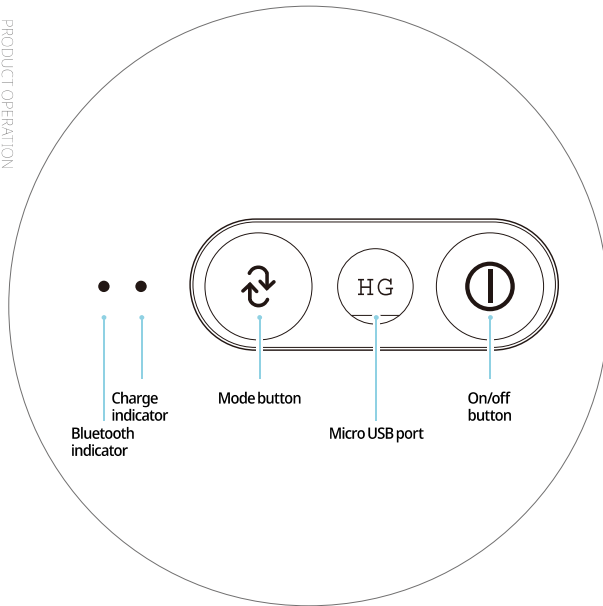EN

Operation

Bluetooth indicator

Connecting or waiting for connection: flash once
 Connected: flash twice

Charge indicator

Low battery: flashing red
 Charging: flashing red
 Fully charged: solid green

Mode button

When the device is turned on, press to switch massaging mode.
 Mode 1
 moving and rotating with heat function on.
 Mode 2
 moving and rotating from front to back with heat function off.
 Mode 3
 moving and rotating from back to front with heat function off.

Micro USB port

Unplug the silicone rubber and connect the cable to the Micro USB port to charge the device. Once fully charged, remove the cable and plug in the rubber.

On/off button

When the device is turned off, press to start up.

When the device is turned on, press to cycle through 3 different levels of intensity (low as default, medium, high).

EN

App connection

※How to download

Scan the QR code to directly download the app
Go to App Market; search "breo" and download

※Operate and use simply according to the guidance given by the app

EN

Charging

- 1 Unplug the silicone rubber.
- 2 Connect the Micro USB to the port on the massager.
- 3 Put the adapter in the socket.

After a "beep", the charge indicator flashes on the control panel.

The indicator lights up continuously when the device is fully charged.

※ The massager and the adapter may become warm during normal charging process.


Time for charging

- 1 Charging takes up to 6 hours.
- 2 Once fully charged, the massager can be used for 10 times (15 minutes a time). (Use time might differ when operating in different temperatures.)
- 3 When the indicator flashes, charge immediately.

Fully Charged

1 Unplug the adapter and remove the Micro USB from the massager.

2 Plug in the silicone rubber.

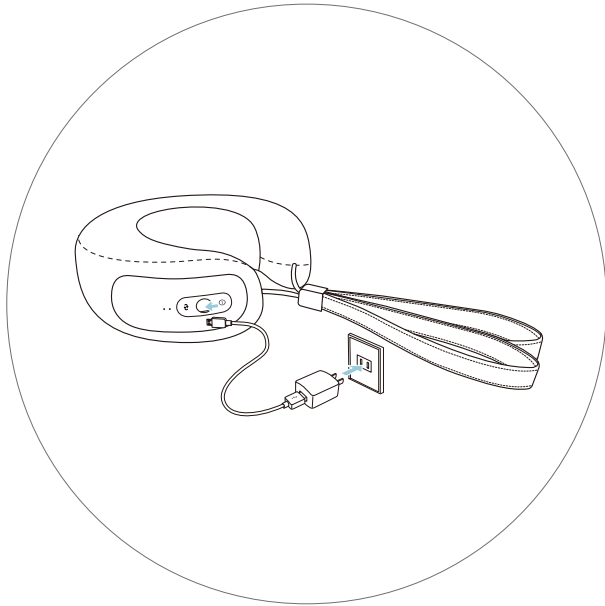17

EN

How to Wear

1 Unlock the magnetic buckle on the bands.

This product fits persons with neck circumference up to 40 cm.


2 Hold both ends of the device and put it on your neck from over to the back of your head.

The control panel might either be on your left or right side.

You can locate different buttons by touching them.

The power button and mode button are both concave. The power button has raised dots on it.

18

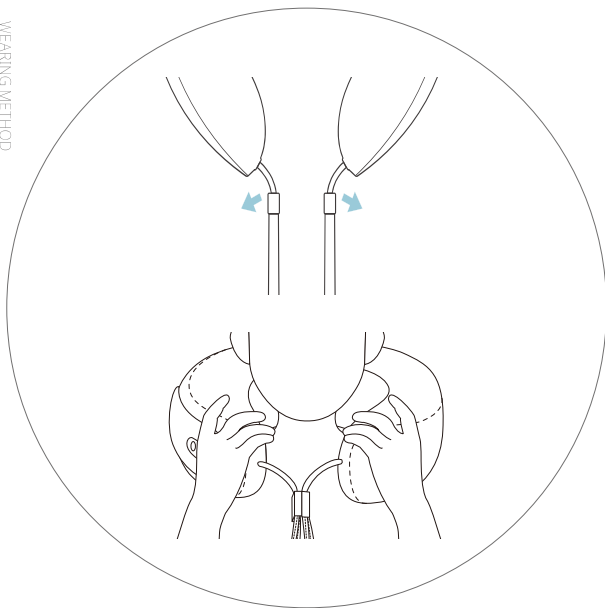**EN****Proper Use**

This device can only be used without cord. Please unplug the adapter after the massager is fully charged.

1 Press the on/off button to switch on.

After a "beep", the control panel will light up and the device will enter Mode 1 as default setting.

2 Adjust the strength of your massage according to your own needs.

Press the on/off button to cycle through low, medium, and high-speed modes.

Press the mode button to choose one of the 3 massaging modes

Mode 1

moving and rotating with heat function on.

Mode 2

moving and rotating from front to back with heat function off.

Mode 3

moving and rotating from back to front with heat function off.

※ **Press and hold the mode button to adjust how the massage head moves.**

※ **You can switch on/off the heat function independently using phone app under all 3 modes.**

3 Wear the massager on your neck. When using this product, please keep your shoulders and neck relaxed, lean back, and relax your whole body. You can pull down the band on both sides to control the pressure.

4 After 10 minutes of use, you will hear a long "beep".
The control panel will light off as the power switch off automatically.

5 Switch off manually

Press and hold the on/off button for more than 2 seconds until you hear a long "beep" to switch off the massager.

6 After use

Unlock the magnetic buckle on the bands. Hold both ends of

the device to take it off.

※ **If you have sensitive skin or this is your first use, observe your neck and skin carefully during use.**

EN

Operation steps for battery removal

- 1 Cut off the external power supply, and remove the product housing;
- 2 Disconnect the power supply to the battery;
- 3 Remove the battery;
- 4 Place the waste battery at a fixed recycling station;
- 5 Please observe local environmental protection laws and waste disposal guidance and protect the environment.

Notes

- Remove the battery out of the appliance before scrapping the appliance;
- The appliance must be disconnected from power supply before removing the battery;
- The battery should be disposed of safely.

Except professionals, no one shall attempt to remove the battery!

EN

Troubleshooting

Please perform the following checks in case of any failure occurring during operation. If the failure cannot be eliminated or the product still cannot work ideally, please cut off power supply and ask a professional technician to repair it. Never attempt to disassemble the product.

Symptom	Possible cause	Basic solution
Power-on failure	Check whether the battery level is too low	If yes, please charge the battery timely
Low massage intensity	Check whether the battery level is low	If yes, please charge the battery timely

If the problem still cannot be solved after troubleshooting, please call our After-sales Service at 400-830-9700

Storage

When the storage period exceeds 3 months, unpack the product, and check and recharge it before use. Do not place the product tightly close to floor, wall and ceiling. The indoor environment should be well-ventilated; keep away from strong sunlight and corrosive gas.

Transport

The product package is provided with simple shockproof measures, therefore is suitable for transport by air, rail, road and sea. Keep away from rain and snow, and avoid inversion and collision.

Check of power supply

Check the power supply; do not power on the product if the supply voltage is beyond the specified range.

Maintenance

Take measures to resist moisture, rats and pollution. Cleaning should be carried out when the product is powered off. After use, put the device in a safe place; do not place any

heavy objects on the device. In case of any failure, please ask a professional technician to repair it. Do not attempt to disassemble the product.

Disposal and Recycle

The product should not be abandoned along with domestic wastes. Dispose of the product according to local environmental laws and guidelines.

Please remember to respect the local regulations: hand in the non-working electrical equipments to an appropriate waste disposal centre. The packaging material is recyclable. Dispose of the packaging in an environmentally friendly manner and make it available for the recyclable material collection-service.

This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

"This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

The statements should be displayed in the user manual: changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or

television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter."

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

This device may not cause interference.

(2) This device must accept any interference, including interference that may cause undesired operation of the device.

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This equipment complies with IC RSS-102 radiation exposure limits set forth for an uncontrolled environment.

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EN

Product name	Neck Massager
Product model	iNeck 3Pro
Rated voltage	5V  1A
Rated power	5W
Battery capacity	2300mAh/1150mAh
Product weight	1000g
Product size	L255×W265×H90(mm)

Shenzhen Breo Technology Co., LTD

19/F, Hisense Southern Building, No.1777 Chuangye Road, Yuehai Street, Nanshan District, Shenzhen 518054, China
support@breocare.com

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Shenzhen Breo Technology Co., Ltd.