

**iSee4**  
EYE MASSAGER

---

USER  
MANUAL

---

**breo.**

---

## **Thank you for purchasing Breo iSee4**

Please read these instructions carefully before using and keep it properly for future reference.

---

## Catalogue


---

Product configuration	3
Function configuration	4
Operation	5
Precautions	7
Troubleshooting and technical performance index	8

## Product configuration

---

The massager includes Master Unit, Power Charger, User Manual .

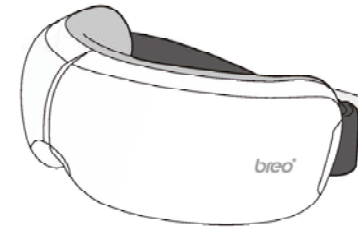

Master Unit



Power Charger




User Manual

## Function configuration

---

Function configuration

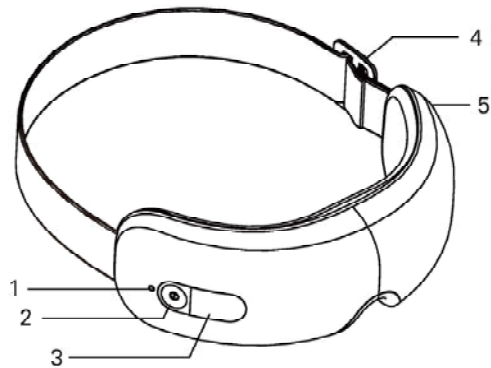


1. Power Supply Indicator
2. On/Off Key
3. Display Area
4. Elastic Adjustment Buckle
5. DC Input Socket

## Operation

---

### Power on

1. Open the lid of the charging port in the bottom of the rear cover of the master unit.
2. Plug the power cord into a wall outlet. We suggest using a surge protected power bar, and plug the other end into the DC jack under the rear cover.
3. The indicator of power on the LCD screen flashes during charging; the flash stops when fully charged.
4. The product cannot be operated while charging (if the product is charged while operating, it will start the self-protection program and stop operating automatically)
5. The time of charging is 2 -3 hours

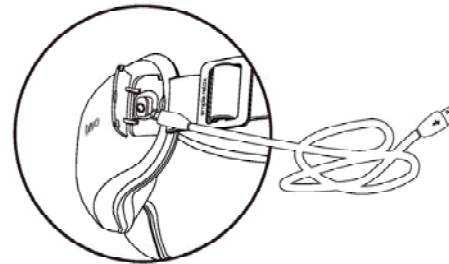


## Operation

---

### Mode Selection

After starting up, the machine's default mode is SLEEP mode; short-press On/Off key to switch next mode.

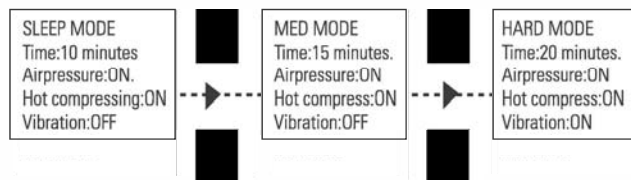

### Shutdown

Long-press On/Off key to shut down the machine, and the machine stops working

## Precautions

---

Please read User Manual carefully to correctly and safely use this product.

1. Please do not use this product while driving.
2. Please do not use this product in the environment with excessively high temperature or humidity.
3. Wear the device and then switch it on to avoid damaging the air bag due to overpressure.
4. Those suffering from eye diseases, such as glaucoma, pearl eye and detachment of retina can not use this product.
5. Those who have had eye surgery can not use this product.
6. The persons suffering from strokes, hypertension, and hypotension should not use this product.
7. Sluggish patients, or people with physical or mental or intellectual disability should not allowed to use this product.
8. Those who have eye surgery can not use this product.
9. Children under 7 years old cannot use this product.
10. Wear the device and then switch it on to avoid damaging the air bag due to overpressure.
11. It is household electric health care product, thus can not be used for diagnosing or curing known disease or other commercial use.
12. If it requires cleaning, wipe it with cotton soaked in alcohol. Do not use water to wash the device.
13. The massager has a hot compressing function, and temperature can reach at 42°C, users not sensitive to change in temperature should be careful in use.
14. Take off your eye contact lenses or glasses before use this product.

## Troubleshooting and technical performance index

---

If there is any problem when you use the device, please follow the instructions below to check. If the problem still exists or the device can't work normally, please turn off the power supply and contact local service center. Do not disassemble or repair the device yourself.

<b>Fault</b>	<b>Possible reason</b>	<b>Solutions</b>
<b>Unable to boot</b>	<b>Whether the cable is in good condition</b>	<b>Check the connection status of the cable</b>
	<b>Whether the battery is low</b>	<b>If yes, charge battery</b>
<b>Air press fomenting is too weak</b>	<b>Whether the battery is low</b>	<b>If yes, Please charge at once</b>

### Technical Performance Index

<b>Power</b>	<b>5V=1A</b>
<b>Power</b>	<b>5W</b>
<b>Ideal environment</b>	<b>Temperature: 5°C-40°C Relative Humidity: ≤80% Atmospheric pressure: 860hpa ~1060hpa</b>
<b>Dimension</b>	<b>Master unit: L 208mm×W 76mm×H 106mm</b>
<b>N.W.</b>	<b>318g</b>